

Full Body Sling



Divided Legs
Reclined Seating
Head Support

Polyester
Model #25010



Polyester Net
Model #25005



Model
#25010



The Full Body Sling is a comfortable sling without reinforcements, handles or padding. It is available in either Polyester or Polyester Net and offers good support for the entire body. The Full Body Sling can be used for most lifting situations and fits most patients.

Versatile Design

The Full Body slings are designed to create a reclined seated position which is extremely helpful if the patient has poor trunk stability and/or a tendency to overstress themselves (spastic performance) or are hypertonic. The reclined position prevents the user from sliding out of the sling. The reinforced head support makes this a suitable sling for lifting patients from a flat surface such as a floor.

Two Fabric Choices

The Full Body slings are available in either Polyester or Polyester Net fabrics. Both can be left behind the patient for a short period of time after the transfer is complete. The low friction Polyester fabric is easy to apply

and remove and the Polyester Net is ideal for most bathing situations.

Easy Care and Easy Fit

Our slings are manufactured to the highest quality standards. They feature reinforced loops and high strength nylon webbing around the sling to make them strong enough to resist even the most active use. Each sling label is colour coded to ensure quick and easy selection of the correct sling and the straps are colour coded for easy determination of required positioning.

HUMAN CARE SWEDEN (HQ)

Årstaängsvägen 21B
117 43 Stockholm
Phone:- +46 8 665 35 00
Fax:- +46 8 665 35 10
info.se@humancaregroup.com

HUMAN CARE CANADA

10-155 Colonnade Road
Ottawa, ON K2E 7K1
Phone: +613.723.6734
Fax: +613.723.1058
info.ca@humancaregroup.com

HUMAN CARE NEDERLAND

Elspeterweg 124
8076 PA Vierhouten
Phone: +31 577 412 171
Fax: +31 577 412 170
info.nl@humancaregroup.com

HUMAN CARE UNITED STATES

8006 Cameron Road, Suite K
Austin, TX 78754
Phone: +512.476.7199
Fax: +512.476.7190
info.us@humancaregroup.com

Accessories



Extension Loops

Used to extend lift straps in order to alter lifting position or to hold leg parts together. Assorted lengths. Part no.:21180.*



Model 21146 Model 25102

Leg Padding and Supports

Ideal for users with sensitive skin or pain. Part no.: 21146-S, M, L Reinforced supports. One size. Part no. 25102



Upholstered Cover

Protects user from hanger bar impact. Not for use in a wet environment. Part no. 50866



Hanger Bars for 2-, 3- or 4-point Suspension

Part no. 57551, 55764 [2- or 3-Point - Quick Release], 55619 [4-Point - Attach to leg hook]

Lifting Situations



Use to/from a bed

Suitable for transfers and repositioning to and from a bed. Can be left behind patient for short period of time.



Use to/from a chair

Suitable for transfers and repositioning to and from a chair. Can be left behind patient for short period of time.



Use to/from the floor


Suitable for transfers from the floor.



Use for showering


Only Polyester Net fabric can be used for this application.

Care and Features

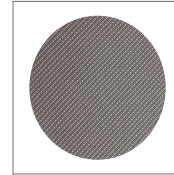


Polyester

Durable, low-friction, easy to clean, apply and remove.

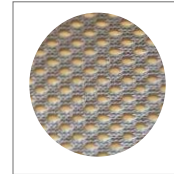

Polyester Net

Soft, thin material that breathes and dries quickly. Ideal for bathing or showering.

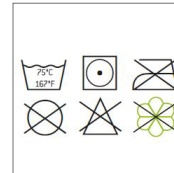

Care

Care instructions located on sling label. Use a hospital grade laundry bag while laundering.



Technical Specifications

Part Number	Size	Colour	Model #25005			Model #25010				
			A	B	C	A	B	C		
Polyester Net	Polyester									
25005-JM	n/a	Junior Medium		22cm / 8.7"	71cm / 28"	80 cm / 31.5"	-	-	-	
25005-S	25010-S	Small		26cm / 10.2"	76cm / 29.9"	86cm / 33.9"	31cm / 12.2"	85cm / 33.5"	90cm / 35.4"	
25005-M	25010-M	Medium		26cm / 10.2"	85cm / 33.5"	95cm / 37.4"	40cm / 15.8"	92cm / 36.2"	91cm / 35.8"	
25005-L	25010-L	Large		34cm / 13.4"	96cm / 37.8"	95cm / 37.4"	40cm / 15.8"	98cm / 38.6"	100cm / 39.4"	
25005-XL	25010-XL	Extra Large		36cm / 14.2"	105cm / 41.3"	106cm / 41.7"	44 cm / 17.3"	105cm / 41.3"	102cm / 40.2"	
n/a	25010-XXL	Double Extra Large		-	-	-	44 cm / 17.3"	111cm / 43.7"	106cm / 41.7"	
n/a	25010-XXXL	Triple Extra Large		-	-	-	44 cm / 17.3"	118cm / 46.5"	108cm / 42.5"	
Maximum User Weight		300kg/661lb								

* Available in assorted lengths: 10cm/3.9" (21180-10), 15cm/5.9" (21180-15), 20cm/7.8" (21180-20), 25cm/9.8" (21180-25), 35cm/13.7" (21180-35)