

Comfort Sling Soft & Comfort Sling High Soft

Sit on Sling

Undivided Legs

Comfortable



Our Comfort Slings are designed to give the patient comfort and a sense of safety. The undivided legs evenly distribute the patient's weight in order to reduce pressure, making this a suitable option for patients with pain. These slings come in two models, one with extra head support.

Unique Fabric

Our Comfort Slings are made of Soft Air Mesh, a material exclusive to Human Care. This soft material promotes less repositioning making it the optimal choice for patients with pain. In addition, the material does not promote skin breakdown and can therefore be left behind the patient for short periods of time.

Functional Design

Our Comfort Slings have no reinforcements or handles and feature undivided legs that evenly distribute the patient's weight. With even distribution, pressure is reduced and patient comfort increased. Our Comfort Slings are available in a medium back and high back

model that offers extra head support. Both slings are considered "sit-on" slings and are applied when the patient is lying down.

Easy Care and Easy Fit

Our slings are manufactured to the highest quality standards. They feature reinforced straps and high strength nylon webbing around the sling to make them strong enough to resist even the most active use. Each sling label is colour coded to ensure quick and easy selection of the correct sling and the straps are colour coded for easy determination of required positioning.



Walking Aids



Healthcare Beds



Lifting Solutions



Hygiene



Convertible Chairs

Accessories



Extension Loops

Used to extend lift straps in order to alter lifting position or to hold leg parts together. Assorted lengths. Part no.: 21180-*



Leg Padding and Supports

Ideal for users with sensitive skin or pain. Part no.: 21146-S, M, L Reinforced supports. One size. Part no. 25102



Upholstered Cover

Protects user from hanger bar impact. Not for use in a wet environment. Part no. 50866



Hanger Bars for 2-, 3- or 4-point Suspension

Part no. 57551, 55764 (2- or 3-Point - Quick Release), 55619 (4-Point - Attach to leg hook)

Lifting Situation



Use to/from a bed

Suitable for transfers and repositioning to and from a bed. Can be left behind patient for short period of time.




Use to/from a chair

Suitable for transfers and repositioning to and from a chair. Can be left behind patient for short period of time.



Use to/from the floor


Suitable for transfers from the floor.

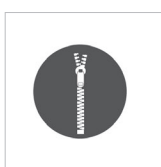

Use for dressing/undressing


Suitable for dressing/undressing applications.

Care and Features

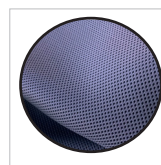

Soft Air Mesh

Lightweight, cool and breathable material. Suitable for sensitive skin.

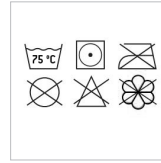

Care

Care instructions located on sling label. Use a hospital grade laundry bag while laundering

Technical specifications

Part Number		Size	Colour	Model #25064			Model #25069			
				A	B	C	A	B	C	
Comfort Sling Soft	Comfort Sing High Soft									Art nr: 25064
25064-JM	25069-JM	Junior Medium	JM	76cm / 30"	73cm / 28.7"	90cm / 35.4"	76cm / 30"	73cm / 28.7"	111cm / 43.7"	
25064-S	25069-S	Small	S	76cm / 30"	80cm / 31.5"	100cm / 39.3"	76cm / 30"	80cm / 31.5"	122cm / 48"	
25064-M	25069-M	Medium	M	86cm / 33.9"	90cm / 35.4"	110cm / 43.3"	86cm / 33.9"	90cm / 35.4"	134cm / 52.8"	
25064-L	25069-L	Large	L	96cm / 37.8"	105cm / 41.3"	116cm / 45.7"	96cm / 37.8"	105cm / 41.3"	139cm / 54.7"	Art nr: 25069
25064-XL	25069-XL	Extra Large	XL	106cm / 41.7"	118cm / 46.5"	120cm / 47.2"	106cm / 41.7"	118cm / 46.5"	145cm / 57.1"	
Maximum User Weight		300kg / 660 lb								

* Available in assorted lengths: 10cm/3.9" (21180-10), 15cm/5.9" (21180-15), 20cm/7.8" (21180-20), 25cm/9.8" (21180-25), 35cm/13.7" (21180-35)



*We improve freedom of movement
and quality of care for people.*

www.humancaregroup.com



HUMAN CARE SWEDEN (HQ)

Årstaängsvägen 21B
117 43 Stockholm
+46 8 665 35 00
info@humancaregroup.com

HUMAN CARE CANADA

10-155 Colonnade Road
ON K2E 7K1 Ottawa
+1 613 723 6734
info.ca@humancaregroup.com

HUMAN CARE THE NETHERLANDS

Elspeterweg 124
8076 PA Vierhouten
+31 577 41 21 71
info.nl@humancaregroup.com

HUMAN CARE UNITED STATES

8006 Cameron rd, Suite K
TX 78754 Austin
+1 512 76 7199
info.us@humancaregroup.com