



Evacuation Sheet

Flame retardant
Low friction
Rescue situation

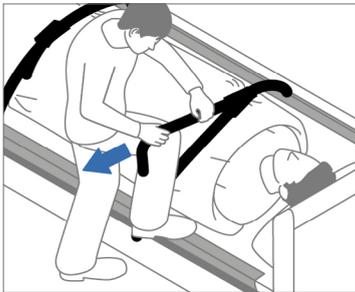

The Evacuation Sheet is used in facilities with clients who are always bedridden. In a rescue situation the Evacuation Sheet can very easily and safely be dragged over floors, tiles, carpeting, down stairs and out of the hazardous area by a single person.

By holding large belt grips at the top or bottom ends and at the centre of the two sides the mattress is pulled out of the bed across the floor into a safe zone. Approved according to DIN1865, Patient Handling equipment for rescue services/ambulances.

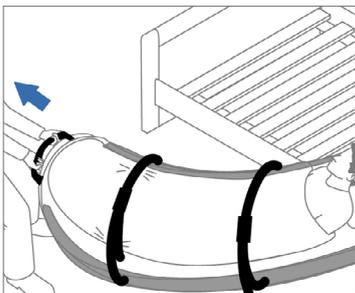


Technical specification

| Article Name | Art. nr. | Dimensions | Material | Max lifting weight |
|------------------|---------------|-----------------------------|---|--------------------|
| Evacuation Sheet | PM-EVA-7053-H | 80 x 200 cm / 34,2" x 78,7" | PVC-coated polyester weave, with ventilation openings, latex-free | 200 kg / 440 lb |



 = 200 kg / 440lb

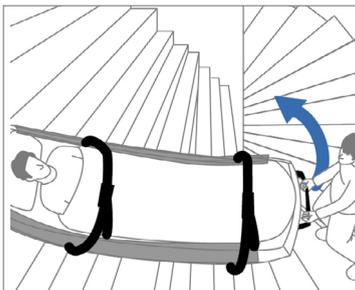

PETERMANN
part of the Human Care Group

These products are distributed by Human Care AB for Petermann GmbH.

Visit www.humancaregroup.com for the latest product information and contact details to national offices and partners.