

Our lift accessories are manufactured to the highest quality standards and are designed to complement and increase the functionality of our lift systems. Simple to install and easy to use, these accessories help to adapt the lift solution to each individual need.

Hanger Bars

We offer multiple options in hanger bars including our newest 2-point quick release series - a lightweight yet durable hanger bar that connects in seconds! They feature an easy-to-read label that indicates the specifications, making it easy to choose the correct hanger bar. With extra large hooks and a high weight capacity, these hanger bars meet most lifting requirements. In addition, we offer 4-point systems to accommodate more robust patients and hanger bars specifically for fixed lift systems.

Medical Grade Scale

Our scale is approved for weighing patients in a medical environment. The scale is quick, accurate and easy to use and compliant with

all relevant standards.

4-Point Powered Spreader Bar

This spreader bar is the optimal choice for hospitals and other health care environments who use slings that are affixed with clips. Controlled by the lift handset, the electric movement is smooth and comfortable for both the client and the caregiver.

Ceiling Lift Trolley

Designed to carry our portable lifts from one room or department to the next, our trolley features a non-slip surface, two directional and two rear locking casters. Easy to assemble it comes complete with a deep basket to carry slings and other lift accessories.

HUMAN CARE SWEDEN (HQ)

Årstaängsvägen 21B
117 43 Stockholm
Phone:- +46 8 665 35 00
Fax:- +46 8 665 35 10
info.se@humancaregroup.com

HUMAN CARE CANADA

10-155 Colonnade Road
Ottawa, ON K2E 7K1
Phone: +613.723.6734
Fax: +613.723.1058
info.ca@humancaregroup.com

HUMAN CARE NEDERLAND

Elspeterweg 124
8076 PA Vierhouten
Phone: +31 577 412 171
Fax: +31 577 412 170
info.nl@humancaregroup.com

HUMAN CARE UNITED STATES

8006 Cameron Road, Suite K
Austin, TX 78754
Phone: +512.476.7199
Fax: +512.476.7190
info.us@humancaregroup.com

Accessories



Hanger Bars

Quick release series.
Altair/Roomer S Part no. 55751, 55764. HeliQ/HeliQ PT Part no. 72754, 72755, 72757



Distance Straps

Used in order to adjust lifting height.



Upholstered Cover

Protects user from hanger bar impact. Hanger bar sold separately. Part no. 50866 & 50866-s



4-Point Suspension

Enables 4-point lifting on HeliQ/HeliQ PT. Part no. 92938



4-Point Suspension

Enables 4-point lifting on Altair/Roomer S. Placed on leg hook. Part no. 55619



Scale

Medical grade scale can be used on all Human Care lifts. Part no. 50590 (Altair/Roomer S), 50591 (HeliQ), 50345 (NA only)



Trolley

Designed to transport Altair/Roomer S. Made of steel (plastic basket). Part no. 70346



4-Point Powered Spreader Bar

For hospitals and other care environments that use clip slings. Three sizes. Part no. 92441, 92442, 92443



Lift Charger

Complete with battery and carry case. For charging lifts. Part no. 71470, 71471, 71472, 71473




Adaptive D-Ring Clip

Allows connection of our portable lifts to existing rail systems. Part no. 53043



Reacher Arm

Assists when hanging/taking down portable lifts so they do not have to be moved by user. Part no. 50494

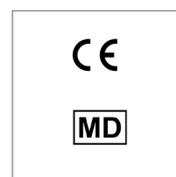

Tested to all relevant standards, including ISO 10535, hoists for the transfer of disabled persons.

Technical Specifications

Hanger Bars						
Part no.	Description	Length	Weight Capacity	Product Weight	Suggested Use	Suspension
55751	Quick Release Series for Altair/ Roomer S	45cm / 17.7"	300kg / 661lb	1.2kg / 2.6lb	Most lifting situations	2, 3 or 4 ¹ -point ¹ with 55619
55764	Quick Release Series for Altair/ Roomer S	55cm / 21.7"	300kg / 661lb	1.5kg / 3.2lb	More robust or sensitive to pressure	2, 3 or 4 ¹ -point ¹ with 55619
72754	Quick Release Series for HeliQ/HeliQ PT	35cm / 13.8"	300kg / 661lb	0.9kg / 2.0lb	Children, Slim Patients	2-point
72755	Quick Release Series for HeliQ/HeliQ PT	45cm / 17.7"	300kg / 661lb	1.1kg / 2.4lb	Most lifting situations	2-point
72757	Quick Release Series for HeliQ/HeliQ PT	60cm / 23.6"	300kg / 661lb	1.3kg / 2.9lb	More robust or sensitive to pressure	2-point
55611	Fixed, for HeliQ/HeliQ PT	35cm / 13.8"	300kg / 661lb	2.0kg / 4.4lb	Children, Slim Patients	2-point
55612	Fixed, for HeliQ/HeliQ PT	55cm / 21.7"	300kg / 661lb	2.0kg / 4.4lb	More robust or sensitive to pressure	2-point
55613	Fixed, for HeliQ/HeliQ PT	60cm / 23.6"	300kg / 661lb	2.0kg / 4.4lb	More robust or sensitive to pressure	2-point
54043	Fixed, for HeliQ/HeliQ PT	45cm / 17.7"	300kg / 661lb	2.0kg / 4.4lb	Most lifting situations	2-point

Hanger Bars for 4-Point Suspension			
Part no.	Dimensions	Weight Capacity	Weight
55619	47cm x 13.4cm / 18.5" x 5.3"	300kg / 661lb	0.8kg / 1.8lb
92938	61cm x 43cm / 24" x 17"	300kg / 661lb	7.2kg / 15.8lb

Accessories			
Part no.	Description	Weight Capacity	Weight
53043	D-Ring	300kg / 661lb	0.1kg / 0.3lb
50494	Reacher Arm	11.2kg/24.6lb	0.5kg / 1.1lb

Continued on next page

Technical Specifications

Trolley	
Part no.	70346
Dimensions (L x W)	63.2cm x 88.4cm /
Height	78.2cm
Height with Handle	88.4cm
Wheel Size Diameter	12.5cm
Basket Size (L x W x H)	47 x 34.5 x 23 cm
Material	Steel / Plastic Basket
Product Weight	14kg
Weight Capacity	20kg / 44.1lb
Compatibility	Altair, Roomer S

Lift Charger	
Part no.	Europe: 71470 Great Britain: 71471 Japan: 71472 North America: 71473
Input	100-240VAC, 1.6A Max, 50 - 60Hz
Output	27.8V, 0.8A, DC
Approvals	IP Class IPX0, CE-marked Compliance to EN-60601-1 (Electrical Safety) and EN-60101-2 (EMC)
Protection	Short circuit, polarization, overheat- ing, Class II against electrical shock
Charging Time	6 hours
Components	Cord, Battery Holder, Battery
Connection	Requires connection to a protective earth power socket

Scale (Global)	
Part no.	50590 (Altair/Roomer S) 50591 (HeliQ)
Dimensions (L x W x D)	12cm x 7cm x 16cm / 4.7" x 2.8" x 6.3"
Display	3cm / 1.5" LCD Panel 5.5 digit display
Battery	6 x 1.5V AAA
Precision Class	2009/23/EC OIML III
Temperature	Ambient: 0°C to +40°C
Approvals	Directive 93/42/EEC, 2009/23/EC CE, FDA D072479/3007074192
Product Weight	1.3kg / 2.9lb
Weight Capacity	300kg / 661.4lb

4-Point Spreader Bar			
	Small	Medium	Large
Part no.	92441 (Manual)	92442 (Electric)	92443 (Electric)
Dimensions (L x W x H)	46.5cm x 45cm x 49.5cm/ 18.3" x 17.7" x 19.4"	46.5cm x 59cm x 49.5cm/ 18.3" x 23.2" x 19.4"	46.5cm x 73cm x 49.5cm/ 18.3" x 28.7" x 19.4"
Input	24V DC, 4A		
Conditions	Ambient Temperature: 10°C to +40°C, Ambient Humidity 30% to 75%, Ambient Pressure 700 hPa to 1060 hPa		
Duty Cycle	-	10%, Maximum 2 minutes ON, 18 minutes OFF	10%, Maximum 2 minutes ON, 18 minutes OFF
Approvals	CE Approved to directive MDD93/42/EEC, In compliance with ISO10535 Models 92442, 92443 IP Rating IPX5		
Product Weight	7.6kg / 16.8lb	10.2kg / 22.5lb	10.6kg / 23.4lb
Weight Capacity	227kg / 500lb		

Scale (North America only)	
Part no.	50348
Dimensions (L x W x D)	16.4cm x 26.5cm x 7cm / 6.5" x 10.4" x 2.7"
Resolution	0 - 200kg = 100g 200 - 320kg = 200g
Battery	4 x LR6/AA
Precision Class	OIML R76 C2/C1.6
Temperature	Ambient: -10°C to +40°C Storage: -40°C to +85°C
Approvals	IP Class: IP54 Directive 90/384/EEC
Product Weight	0.7kg / 1.5lb
Weight Capacity	320kg / 705lb

Upholstered Cover			
Part no.	50866-s 50866	Material	Polyester
Dimensions (L x W)	26cm x 23cm / 10.2" x 9.1" 36cm x 23cm / 14.2" x 9.1"	Hanger Bar Compatibility	55673 55751, 55764, 72755, 72757
Warnings	Not for use in a wet environment		